




Product - Information

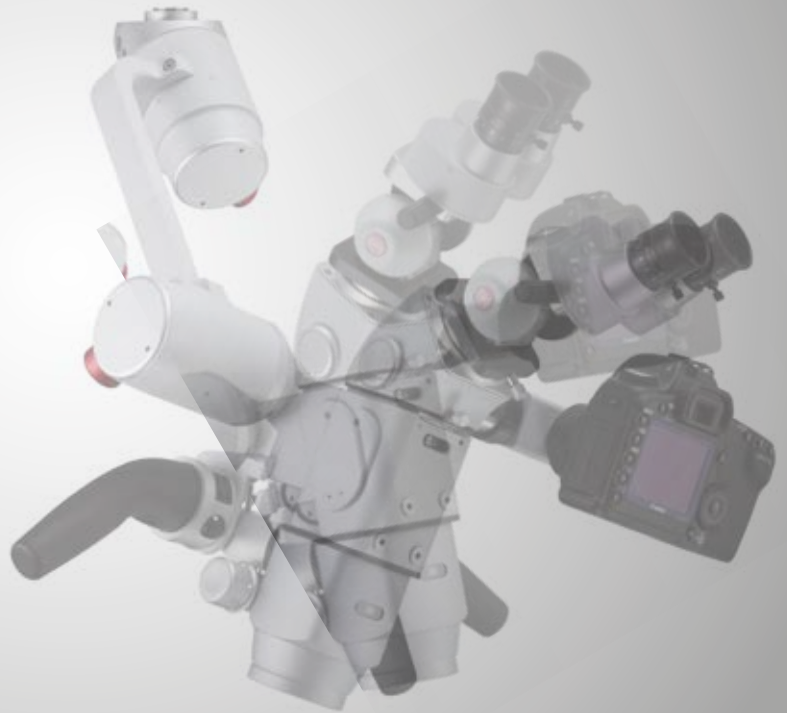

Kaps - 45° Ergonomic Wedge with Beam Splitter and Binocular Turntable

An unbeatable combination for a high level of convenience at the workplace.

Dentists, in particular, often complain about back or neck problems. Thanks to the Kaps 45° ergonomic wedge, including binocular turntable, these are now a thing of the past.


The user always sits comfortably in an upright, ergonomically correct position and can work in a much more relaxed way.

Even with varying working positions, or if the microscope is pivoted away to the side, the user always keeps the patient in focus. A simple hand movement and the preferred working position can be assumed again. The neck and back are noticeably relieved.



45° Ergonomics Handle with Beam Splitter

The 45° ergonomics handle with beam splitter is simply inserted between the magnifying changer/zoom and binocular turntable.

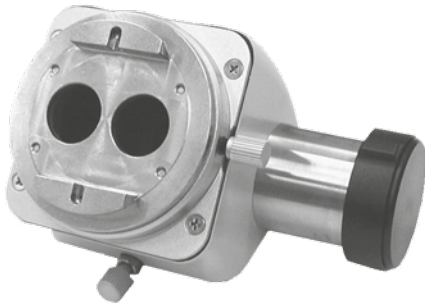

45° ergonomics handle unilateral 50:50

➔ Order no 101-990-100

45° ergonomics handle unilateral 80:20

➔ Order no 101-991-100



45° ergonomics handle bilateral 50:50 / 80:20

➔ Order no 101-993-100

Binocular Turntable

- The low height of just 20 mm ensures a minimal extension of the microscope head
- It weighs only 240 g. This means that no harmful interference occurs during handling of the microscope
- The rotation angle is $\pm 24^\circ$ and is rasterized in $4 \times 6^\circ$ steps
- With the **clamping screw**, locking of the position is possible in any position
- It is always possible to adapt all Kaps microscopes at a later point
- An adaptation of competitor models is also possible

➔ Order no. 101.989-000

